

DC Dishcharge Experiment Lab Report

Khush Thakor, Chad Anglemyer, Tristan Lewis, Matthew Ryan

11/19/2022

Purpose

Safely ionize air into a plasma and investigate the relationship between breakdown voltage and pressure multiplied by distance between cathode and anode.

Abstract

What is plasma? Plasma is the ionization of a gas, meaning we are separating the protons and electrons from one another. Plasma often produces a visible light from electrons accelerating.

Figure 1 is an image of our experiment layout. Additionally, we do have a vacuum pump located on the floor to the left of the pressure gauge. The way our experiment works is by using our vacuum pump to extract air from our gas tube which results in a lower pressure in our glass tube. We installed a fine adjust valve (pressure valve 2 in the figure) to keep the air

between pressure valve 1 and pressure 2 as control air which we could have enter our glass tube to maintain a stable pressure. However, were not able to adjust the valve efficiently to let minimal air into the tube in order to keep stable pressure. After the vacuum tube was at our ideal pressure, we then turned on the high voltage power supply and carefully adding voltage until we saw a breakdown of plasma in the vacuum tube.

During preliminary investigation we had a break down of plasma at 23.35 cm distance between cathode and anode approximately 600 millitorr pressure and 1.02 kilovolts of power from anode and cathode. With our given set up we are able to change 3 variables of our experiment, the distance between our anode and cathode, pressure in our glass tube and voltage applied to anode and cathode. Initial research prompted experimenting with the idea of Paschen curve. According to "Foundations of dc plasma sources" by J.T Gud-

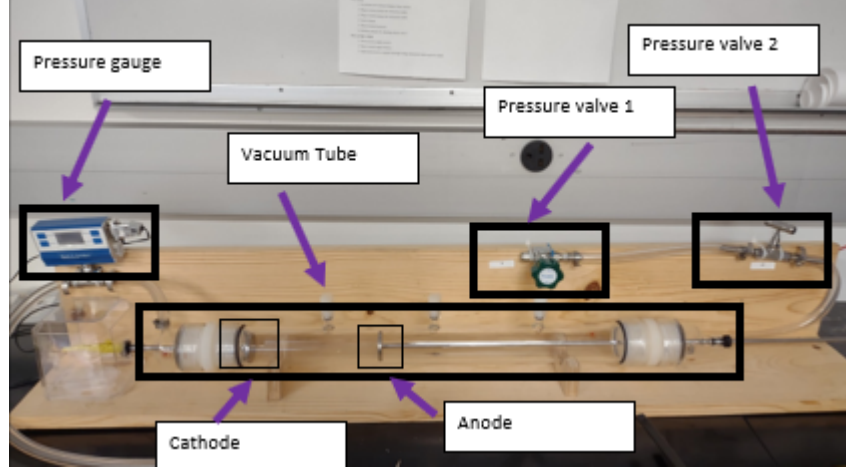


Figure 1: DC Plasma Formation Experiment Diagram

mundsson and A. Hecimovic, “This curve shows that for a fixed discharge length L there is an optimum pressure for plasma break-down.” The source also states, “At lower pressures the ionization process is ineffective due to the low electron -neutral collision probability while at higher pressures elastic collision prevent the electrons from reaching high enough energy.” Explaining how at lower pressure the distance between cathode and anode needs to be longer in order to produce a discharge similar to that of high pressure and low distance between cathode and anode.

Given this correlation we had a working hypothesis relating all com-

ponents distance, pressure, voltage. We were able to say that pressure multiplied by distance is related to voltage, by decreasing our pressure-distance ratio, leads to less voltage required to produce plasma. Once we had this working hypothesis, we plotted a total of 7 distances and 3 data points at all distances besides the lowest distance of 1.01cm we acquired 5 points. After finding all this data we were able to relate the pressure multiplied by distance to voltage by the given equation. $y = -1E - 09x^3 + 3E - 06Ex^2 - .0013x + 0.8273$, where Y is the voltage and X is pressure multiplied by distance.

Procedures

1. Remove clutter around and close to the plasma tube and equipment
2. Check connections on the plasma tube and the power supply for no loose wires or other objects
3. Measure distance between electrodes for the tube and adjust to 23.5 cm
4. Check viewing glass in oil pump and make sure oil level is $\frac{3}{4}$ full. If oil is needed follow instructions in manual to refuel pump
5. Plug in cable to the power supply and then into the wall socket
6. Plug in the pump and pressure gauge
7. Switch on the motor pump
8. Wait for pressure to reach 350 mT
9. After pressure is reached, turn on switch 1 on the front of the power supply, then switch 2 on the front of the power supply
10. Turn voltage knob slowly, making note per interval of 0.1 volts
11. Visually watch the glass tube for signs of plasma, once initial plasma is shown record pressure reading on pressure gauge and voltage reading on power supply
12. Slowly turn voltage knob on power supply slowly down towards zero.
13. Turn off switch 2 on front power supply and then switch 1 on front power supply
14. Release valve 1 and slowly release valve 2 and make sure that tube pressure reaches 77- torr to prevent oil pump air from entering experiment,
15. Shut off pump
16. Unplug power supply and pump

17. Repeat steps from step 1 for 3 times for each electrode distance and then change the electrode distance to 18.5 cm, 13.51 cm, 8.52 cm, 3.5 cm and finally 1 cm

Data

| Electrode Distance (cm) | Pressure Reading (mTorr) | Plasma Breakdown Voltage (kv) |
|-------------------------|--------------------------|-------------------------------|
| 23.50 | 365 | 1.02 |
| 23.50 | 330 | 1.07 |
| 23.50 | 315 | 1.07 |
| 28.50 | 345 | 1.11 |
| 28.50 | 335 | 1.15 |
| 28.50 | 341 | 1.14 |
| 18.52 | 345 | 1.02 |
| 18.52 | 345 | 0.96 |
| 18.52 | 332 | 0.93 |
| 13.51 | 370 | 0.88 |
| 13.51 | 337 | 0.68 |
| 13.51 | 327 | 0.7 |
| 8.52 | 355 | 0.85 |
| 8.52 | 333 | 0.63 |
| 8.52 | 306 | 0.66 |
| 3.51 | 321 | 0.6 |
| 3.51 | 260 | 0.65 |
| 3.51 | 255 | 0.65 |
| 1.01 | 320 | 0.6 |
| 1.01 | 281 | 0.77 |
| 1.01 | 768 | 0.38 |
| 1.01 | 317 | 1.27 |
| 1.01 | 273 | 0.61 |

Calculations

$$1\text{Pa} = 0.33\text{mT}$$

| Electrode Distance (cm) | Pressure Reading (mTorr) | Plasma Breakdown Voltage (kv) | P (Pa) | $P \cdot D$ (mTorr*cm) | $P \cdot D$ (Pa*cm) |
|-------------------------|--------------------------|-------------------------------|--------|------------------------|---------------------|
| 23.5 | 365 | 1.02 | 48.55 | 8578 | 1140 |
| 23.5 | 330 | 1.07 | 43.89 | 7755 | 1030 |
| 23.5 | 315 | 1.07 | 41.9 | 7403 | 985 |
| 28.5 | 345 | 1.11 | 45.89 | 9833 | 1310 |
| 28.5 | 335 | 1.15 | 44.56 | 9548 | 1270 |
| 28.5 | 341 | 1.14 | 45.35 | 9719 | 1290 |
| 18.52 | 345 | 1.02 | 45.89 | 6389 | 850 |
| 18.52 | 345 | 0.96 | 45.89 | 6389 | 850 |
| 18.52 | 332 | 0.93 | 44.16 | 6149 | 818 |
| 13.51 | 370 | 0.88 | 49.21 | 4999 | 665 |
| 13.51 | 337 | 0.68 | 44.82 | 4553 | 606 |
| 13.51 | 327 | 0.7 | 43.49 | 4418 | 588 |
| 8.52 | 355 | 0.85 | 47.22 | 3025 | 402 |
| 8.52 | 333 | 0.63 | 44.29 | 2837 | 377 |
| 8.52 | 306 | 0.66 | 40.7 | 2607 | 347 |
| 3.51 | 321 | 0.6 | 42.69 | 1127 | 150 |
| 3.51 | 260 | 0.65 | 34.58 | 912.6 | 121 |
| 3.51 | 255 | 0.65 | 33.92 | 895.1 | 119 |
| 1.01 | 320 | 0.6 | 42.56 | 323.2 | 43 |
| 1.01 | 281 | 0.77 | 37.37 | 283.8 | 37.7 |
| 1.01 | 317 | 1.27 | 42.16 | 320.2 | 42.6 |
| 1.01 | 273 | 0.61 | 36.31 | 275.7 | 36.7 |

Error Discussion

Uncertainty and Error

The tools we used for measurements, a ruler down to the mm allowing for a ± 1 mm error. We also used the Pressure gage which to the millitorr allowing for a ± 1 millitorr. The object with the most error was the power supply which a ± 0.01 kv error, shown in Figure 2.

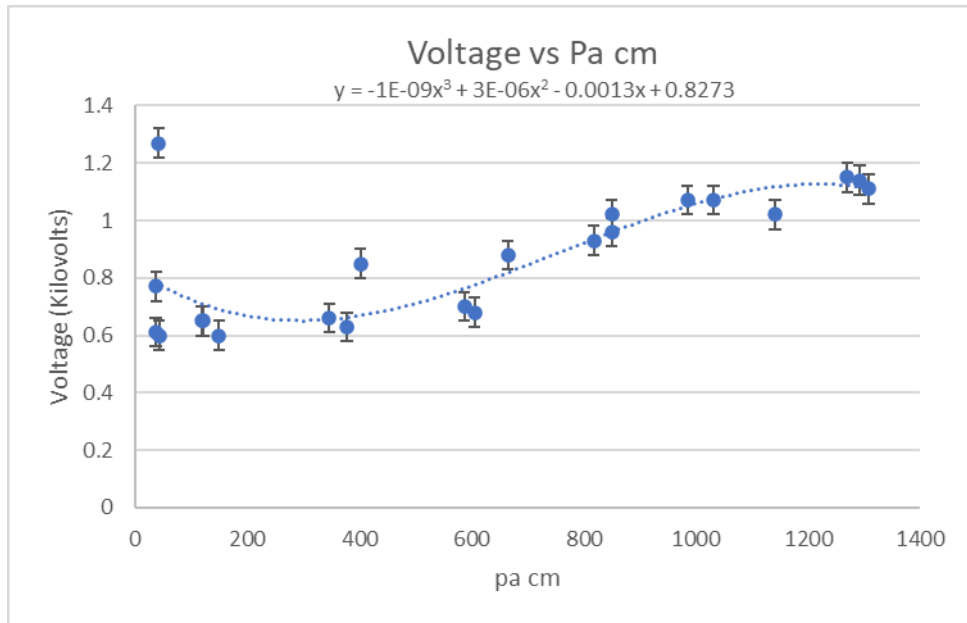


Figure 2: Voltage over Pa cm Graph

Mistakes/Blunders

In our experiment we are not able to keep the pressure constant in our plasma tube which lead to vary in our data however did create plot our data with pressure multiplied by distance.

Result

During this experiment it was quite apparent that we had close to zero control over a stable pressure. Another topic of reflection would be to organize a another experiment with more control on our variables in order to which shows the Pachens Curve. With our data collected from the experiment we found the relationship between pressure, electrode distance and voltage. The relationship found was that as the pressure-distance ratio decreased, the lower the voltage was required to make plasma. This matched our hypothesis, and the experiment data supported it.

Result Statement

The calculations show that when the electrode distance decreased, and the pressure decreased in the containment tube, the required voltage to obtain plasma also decreased.

Future Experiments and Possible Alterations

Looking into the possibilities of the experiment, one of the biggest prospects is a different vacuum tube setup. Possible shapes for a new containment chamber can be dome like or if the tube shape is retained, possible concave or convex portions of the tube can alter the experiment. To further the observation ability of the plasma different apparatuses can be placed into the containment chamber for voltage application or magnetic application. An example of this can be a coil to change the movement of the plasma, and to further observe

the effect of magnetism on the plasma you can input magnets into or around the containment chamber. With the current setup it is possible to analyze the same pressure-distance ratio and how that affects plasma, however now we can run the experiment with different gases, for example, the element Argon. This will allow for a new range of data to create and compare it to the other tests. Looking into the issues we had, something that needs to be resolved is the lack of pressure control. A possible solution to this would be access to a pressure valve between the vacuum pump and the containment tube. This would allow the user of the system to isolate the system of the plasma experiment and hopefully have a stable pressure inside the system. The only concern with this is that the new stable vacuum would have to be carefully released when the desired tests have been completed